

 **QLICKSMART®**  
**Snapit**  
Multi-Use | Ampoule Opener



Work Safe

Play Safe Tomorrow.

# ▲ The underestimated cut injury from glass ampoule

One in three healthcare workers has experienced a cut when opening a glass ampoule by hand. Regrettably, many users and clinical administrators have underestimated the likelihood of these injuries.

Cuts from glass ampoules can happen to anyone, including you, no matter how experienced you are in opening ampoules. Even a minor ampoule cut can restrict your clinical activities. More serious injuries such as tendon lacerations have been reported.

Repeated personal injuries whilst opening ampoules, led an Australian registered nurse to invent the award-winning Qlicksmart SnapIT ampoule opener.


## ▲ The award-winning ampoule opener



Qlicksmart's SnapIT is designed to protect users from injury when an ampoule exposes a sharp edge or shatters while being opened.



The ampoule openers are simple to use with a quick "insert and snap" action.

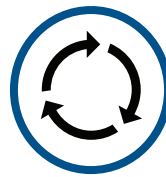

Being reusable, SnapIT is a cost-effective solution to open glass ampoules.



The ampoule opener has won multiple awards, which demonstrate its practicality and effectiveness in reducing glass ampoule injuries.

### Lite Edition



- SnapIT Lite is perfect for hanging on a lanyard or carrying in a pocket for easy, on-the-go use.
- The ampoule opener is made from a sturdy and high quality plastic.
- A hexagonal anti-roll collar ensures the SnapIT Lite ampoule opener remains stable when placed on a table or trolley

### Personal Edition



- Made from high quality aluminum with two colour choices.
- The Personal edition is perfect for hanging on a lanyard or carrying in a pocket for easy, on-the-go use.
- Autoclavable if required\*.

### Trolley Edition



- Made from high quality aluminum.
- The Trolley edition is available in three sizes (Regular, Large, and Extra Large); all with a flat side that stops it from rolling and potentially falling off a table or trolley.
- Autoclavable if required\*.

Product edition	Personal	Trolley			Lite	
Product code	SN-01R	TE-01R	TE-01L	TE-01XL	SN-02R	SN-02L
Product size	Regular	Regular	Large	Extra Large	Regular	Large
Compatible ampoule sizes*	1-2ml, 5-10ml, 10-15ml	1-2ml, 5-10ml, 10-15ml	5-10ml, 10-15ml, 20-25ml	20-25ml, 25-30ml	1-2ml, 5-10ml, 10-15ml	5-10ml, 10-15ml, 20-25ml
Size per unit	Length: 80mm Width: 17mm	Length: 80mm Width: 25mm	Length: 80mm Width: 29mm	Length: 80mm Width: 25mm	Length: 80mm Width: 13mm	Length: 80mm Width: 18mm
Weight per unit	34g	50g	74g	66g	11g	13g
Colour range	Blue, Purple	Black	Red	Blue	Metallic Grey, Purple, Teal, Pink, Blue.	Burgundy
GMDN code	10098					
Material	Anodised aluminium (body), Silicon (O-ring), Stainless steel (keyring and spring).				Plastic (body), Silicon (O-ring), Stainless steel (keyring and spring).	
Latex	No					
DEHP	No					
Sterile	No, but autoclavable to 130 °				Not autoclavable.	
Re-usable	Multi-use. O-ring can be replaced if lost or broken.					
Regulatory	USA FDA listed					

\* Compatibility depends on the ampoule's geometry and the diameter of the ampoule lid.  
Ampoule of the same size can have a different geometry and ampoule lid diameter. The ampoule size in ml is only used as a guide.

## How to use SnapIT

**SnapIT**  
Multi-Use | Ampoule Opener




1. Align the dots if present. Gently insert all of the ampoule lid into the SnapIT Lite. The ampoule lid becomes locked in at the neck.



2. With little effort, lever away from the dot to snap off the lid. The lid should remain held in the SnapIT.

**NOTE:** If your standard practice requires you to turn the dots away from your body, make sure that both dots are still aligned after turning and before snapping away from dot.

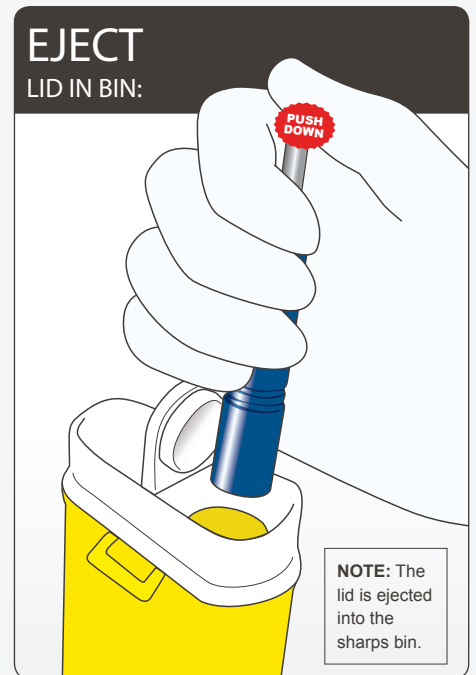

3. Carefully aim over a sharps bin and eject the lid by pushing down on the rod with your thumb.

**NOTE:** The lid is ejected into the sharps bin.



### CLEANING YOUR SNAPIT:

4. Clean your SnapIT regularly. Remove O-ring as shown to disassemble the parts and clean with water or hospital grade cleaner. Do this in a clean dish to avoid loss of parts. Dry completely then re-assemble.

**NOTE:** If the ampoule lid breaks, tap the SnapIT over a sharps bin to remove excess glass pieces then rinse under running water before disassembling.

#### SIZES:

If the ampoule lid is too big, it will not fit through the entrance. Use a larger SnapIT.  
If the ampoule is too small, the SnapIT will not hold the top securely. Use a smaller SnapIT.

## Importer / Distributor :



MEDZORG Nederland BV  
Satellietbaan 14 B  
2181 MH Hillegom NL



0252 - 536200



orders@medzorg.nl  
info@qlicksmart.nl



[www.qlicksmart.com](http://www.qlicksmart.com)

## Sub-Distributor :

